

INDICATIONS

IMPORTANT: Begin by brushing and rinsing your teeth thoroughly. All food debris must be removed from brackets before using your OraKlenz Braces Care System.

Your dentist will demonstrate how each tray is inserted and removed. Trays are separated into quadrant sections (4 trays total). Each quadrant section must be inserted separately.

GEL APPLICATION:

OraKlenz Trays are designed to hold the gel right where it is needed. This ensures that the gel penetrates below the tissue and into the pockets. Place the trays open side up and apply a line of OraKlenz Gel across all the tooth indentations. **DO NOT** over fill the tray with gel. A small amount in each tooth indentation is all that is needed.

NOTE: To avoid contamination, **DO NOT** touch the sides of the tray with the applicator.



TRAY INSERTION:


To insert tray, approach from the bracket side of teeth at a 45 degree angle. Position the tray so that the tray edge bypasses the bracket, then roll tray up and over the bracket. Always check to make sure the trays are seated firmly in place.

CAUTION: Insertion of upper and lower trays should be done using **FINGER PRESSURE ONLY** – thumbs to seat the upper and index finger to seat the lower. **Trays should NEVER** be “bitten” into place.



TRAY REMOVAL:

To remove tray, press finger against tongue side, rolling it up, over and away from the brackets. **NOTE:** Take care to ensure the tray does not “grab” the bracket as it is removed.

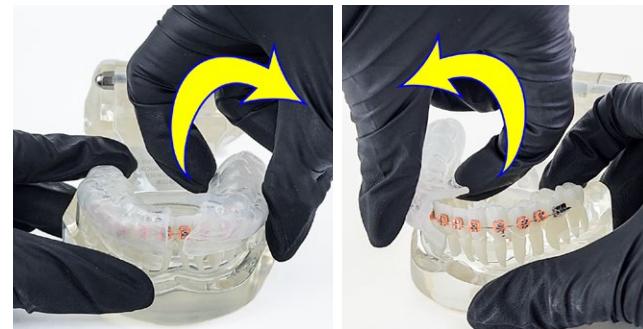

DAILY TREATMENT:

Use your thumb and index finger to rub and press the outer surface of the trays. Then, gently bite down so that the upper and lower trays make contact.

Keep the trays in place for 15-20 minutes. Depending on the severity of your periodontal condition, this process should be repeated 1-3 times a day. NOTE: We have found that the best time to perform your daily treatments is after mealtimes, shower/bath times, TV relaxation time, or just before bedtime.

AFTER TREATMENT:

When treatment is completed, remove the trays and spit any excess gel into the sink. Rinse the trays under cool water – **NEVER USE HOT WATER** -- and brush each tray with a soft toothbrush. Allow trays to dry until the next treatment.

At least once a week, trays should be cleaned with an appliance cleaner, e.g. EverSmile White Foam Appliance Cleaner. When properly used, 1 tube of gel should last about 4 weeks. Always keep the gel in a cool dark place, such as a medicine cabinet or drawer. **DO NOT EXPOSE TO LIGHT.**

Doctor's Instructions:

Wear Time – 15 to 20 minutes

☐ 3 Times Daily ☐ 2 Times Daily ☐ Once per Day

Special Instructions:

Dr. _____ Phone _____

Important Note: OraKlenz does not replace regular brushing and flossing – If you feel any irregular or adverse reaction contact your dentist immediately.

FABRICATED BY:



SML® - Space Maintainers Laboratories
9129 Lurline Ave, Chatsworth, CA 91311
1-800-432-3270
www.SMLglobal.com



Patient Instructions

The OraKlenz™ Braces Care System

Congratulations! You've taken a giant step towards improving your dental health. Your dentist has prescribed the OraKlenz Braces Care System. These trays deliver healing OraKlenz Gel directly to the source of your infection. Please follow your dentist's treatment instructions carefully.



CAUTION: Federal law requires a written prescription, signed by a licensed medical professional.



REV.C.MAY.2018